INSTRUCTIONS FOR USE

185 044 542 Rev B

CONTENTS

English

Français (French)


Espanol (Spanish)

FISHER & PAYKEL HEALTHCARE

INTERNATIONAL	PO Box 14 348, Papanui, Auckland 1741, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fpicare.com Website: www.fpicare.com
USA	Tel: 1800 446 3908 or +1 949 470 3900 Fax: +1 949 470 3933
AUSTRALIA	Tel: +61 3 9879 5022 Fax: +61 3 9879 5232
AUSTRIA	Tel: 0800 29 31 23 Fax: 0800 29 31 22
CHINA	Tel: +86 (20) 3205 3486 Fax: +86 (20) 3205 3495
FRANCE	Tel: +33 (0)1 64 46 52 01 Fax: +33 (0)1 64 46 52 21
GERMANY	Tel: +49 7182 93777 0 Fax: +49 7182 93777 99
INDIA	Tel: +91 80 4123 6041 Fax: +91 80 4123 6044
IRISH REPUBLIC	Tel: +1800 409 011
ITALY	Tel: +39 06 7839 2939 Fax: +39 06 7814 7709
JAPAN	Tel: +81 (3) 3661 7205 Fax: +81 (3) 3661 7206
NETHERLANDS	Tel: +31 78 644 0924 Fax: +31 78 644 0914
SPAIN	Tel: +34 902 013 346 Fax: +34 902 013 379
SWITZERLAND	Tel: 0800 83 47 63 Fax: 0800 83 47 54
UK (EU Authorized Representative)	Tel: +44 1628 626 136 Fax: +44 1628 626 146
NTHN IRELAND	Tel: +0800 132 189

FlexiFit™ 432 Full Face Mask Instructions

The FlexiFit™ 432 Full Face Mask is intended for multiple patient or single patient adult use by individuals who have been diagnosed by a physician as requiring CPAP or Bi-level ventilator treatment in the home, hospital or other clinical setting.

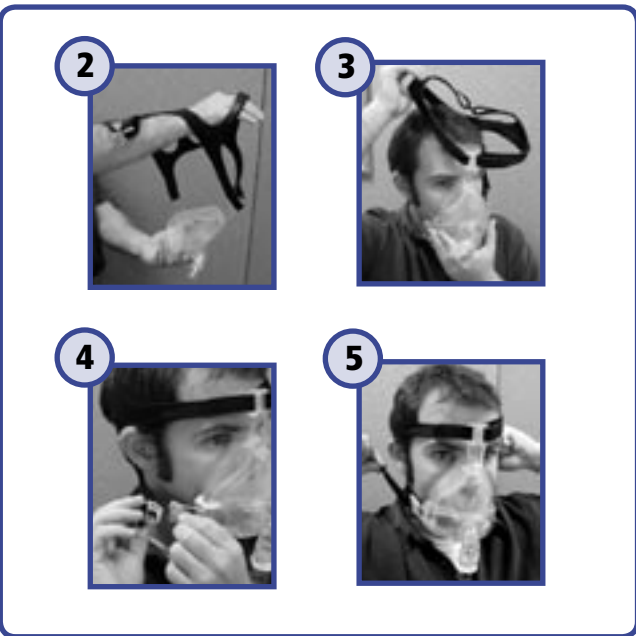FITTING YOUR MASK

Refer to the mask assembly diagram

- Remove plastic cover from Silicone Seal.
- Unclip the Headgear (K) from the Glider™ strap (E) by pulling the Release Cord (I) down. Holding your Mask Base (A) in one hand, spread the Headgear with your free hand just enough to fit comfortably over your head.
- Lightly hold the mask against your face, with the bottom of the seal under your chin, and spread the Headgear over your head with your free hand.
- Re-attach the Headgear Clip (H) to the Glider™ strap.
- Gently adjust the Top Horizontal straps, then the Lower Horizontal straps, then the Top Crown straps. Ensure the top of the seal remains centered on the bridge of the nose, and the bottom of the seal remains fixed under the chin. Ensure that the Glider™ strap is centered.
- Connect the CPAP machine hose running at the prescribed pressure to the mask.

FITTING TIPS

- If necessary pull the mask slightly away from your face while the straps are still attached. This allows the silicone seal to inflate. Put the mask back on your face.
- If leaks occur around the upper half of the mask gently tighten the top horizontal straps. It is not necessary for the Forehead Cushion to touch your forehead.
- If leaks occur around the lower half of the mask gently tighten the lower horizontal straps. (Avoid over tightening the straps.)
- When you need to remove your mask unclip the Headgear from the Glider™ strap by pulling the Release Cord forward or down.

CLEANING YOUR MASK

To ensure continuing safety and comfort for your FlexiFit™ 432 Full Face Mask:

Before first use:

- Hand wash the Headgear in pure soap dissolved in lukewarm water. Do not soak for longer than 10 minutes, then rinse in fresh water and allow to air dry out of direct sunlight. It is not necessary to remove the Clip (H) from the headgear when washing.

Before each use:

- Inspect the mask for deterioration. Do not use if aged or damaged.
- Inspect the Non-Rebreathing Valve (J). Ensure the Silicone flaps are down when the CPAP is turned off.


After each use:

- Disassemble the Silicone Seal (C), Foam Cushion (D), Elbow (B), Non-Rebreathing Valve (J), and the Headgear (K) from the Mask Base (A).
- Wash the Silicone Seal, Elbow, Valve and Mask Base in pure soap dissolved in lukewarm water. Do not soak for longer than 10 minutes, then rinse in fresh water.
- Wipe the Foam Cushion (D) with a damp cloth (do not immerse the Foam Cushion in water).
- Leave all parts to dry, away from direct sun light before reassembling.

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- Caution:**
- Do not soak the mask for longer than 10 minutes.
 - Do not wash the mask in the dishwasher.
 - Do not clean the mask with products containing alcohol, chlorine, antibacterial agents or bleach.
 - Do not store the mask in direct sunlight.
 - The above actions may deteriorate or damage the mask and shorten its life. If your FlexiFit™ 432 Full Face Mask weakens or cracks, replace immediately.

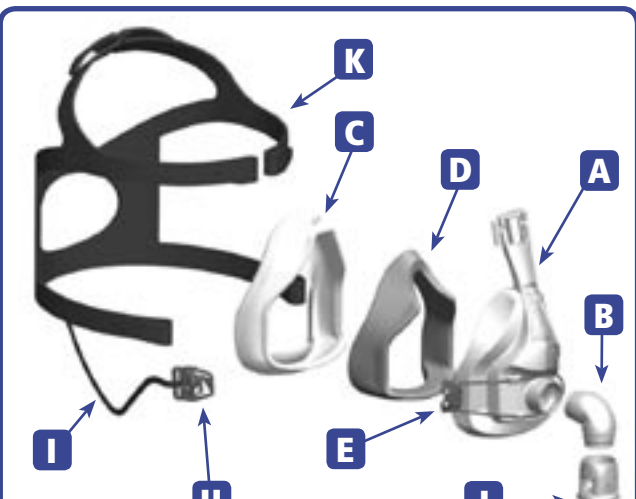
ASSEMBLING YOUR MASK

Your FlexiFit™ 432 Full Face Mask is fully assembled and ready to use. However, when you need to clean your mask follow these instructions for reassembly.

- Place Silicone Seal (C) face down on a flat surface. With the thick side of the Foam Cushion (D) facing down, place onto Silicone Seal.
- Place Silicone Seal (C) and Foam Cushion (D) into the Mask Base (A).
- Push the Non-Rebreathing Valve (J) onto the Elbow (B), and attach to Mask Base.
- Lay the Headgear (K) flat with the label side down. Position the assembled mask over the Headgear and attach the four straps in to the corresponding slots in the Mask Base and the Glider™ strap (E). This can be done without undoing the Velcro™ Tabs by sliding the Headgear into the slots.

The FlexiFit™ 432 Full Face Mask and accessories are Latex free. Many people find CPAP therapy causes dryness of the nose and throat. Heated humidification can substantially reduce these symptoms. Contact your supplier and ask about ThermoSmart™.

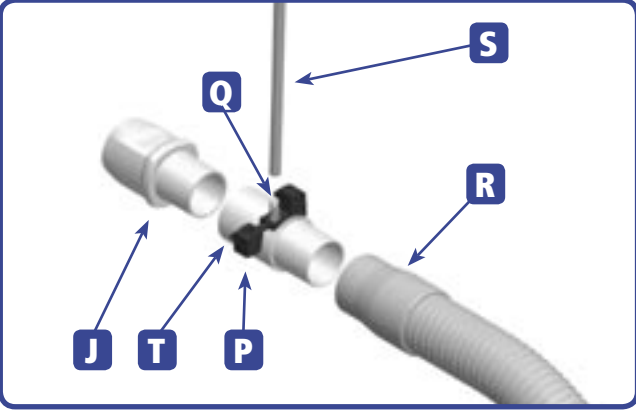OXYGEN / PRESSURE PORT CONNECTOR

If pressure readings and/or additional oxygen are required, an Oxygen/Pressure Port Connector is available. (Code 900HC452)

- Connect the large end of the Connector (T) to the Non-Rebreathing Valve (J) on the mask and the small end of the Connector to the Breathing Tube (R).
- Open the Port Caps (P) and place the Pressure and/or Oxygen Tube (S) firmly over the Port(s) (Q).

The two Ports are identical; either Port may be used for oxygen or pressure measurement. Only remove the Port Cap from Ports that are being used.

- Note: At a fixed flow rate of supplemental oxygen flow, the inhaled oxygen concentration will vary, depending on the pressure settings, patient breathing pattern, mask selection, and the leak rate.
- Warning: If oxygen is used with this CPAP machine, the oxygen flow must be turned off when the CPAP machine is not operating.
- Explanation of Warning: When the CPAP device is not in operation, and the oxygen flow is left on, oxygen delivered into the ventilator tubing may accumulate within the CPAP machine enclosure. Oxygen accumulated in the CPAP machine enclosure will create a risk of fire.
- Warning: Do not smoke while using or in the vicinity of someone using oxygen with this mask.
- Warning: Do not use Polyvinyl Chloride (PVC) Pressure and/or Oxygen tubing. Use silicone tubing instead.
- Explanation of Warning: Contact with PVC can cause premature cracking of ports.

WARNINGS

- This Mask should be used only with CPAP or Bi-Level equipment recommended by your physician or respiratory therapist. A mask should not be used unless the CPAP or Bi-Level machine is turned on and operating properly. The exhalation ports associated with the mask should never be blocked.
- Explanation of Warning: CPAP and Bi-Level Machines are intended to be used with special masks which have exhalation ports to allow continuous flow of air out of the mask. When the CPAP or Bi-Level machine is turned on and functioning properly, new air from the CPAP or Bi-Level machine flushes the exhaled air out through the exhalation ports. However, when the CPAP or Bi-Level machine is not operating, enough fresh air will not be provided through the mask, and exhaled air may be re-breathed. Rebreathing of exhaled air for longer than several minutes can in some circumstances lead to suffocations.
- At low CPAP pressures the flow through the exhalation ports may be inadequate to clear all exhaled gas from the tubing. Some rebreathing may occur.
- Operating pressure range is 3-25 cm H₂O.
- Valve activation/deactivation pressure is 2 cm H₂O.
- Mask deadspace is approximately 115 cc including the swivel and elbow. This volume will vary slightly depending on individual facial features.
- Do not use the mask without the Non-Rebreathing Valve in place. This valve ensures CO₂ rebreathing is kept to a minimum in the event of power failure.
- Do not use the mask if vomiting or feeling nauseous.
- Discontinue using the mask if there is an allergic reaction to any part of the mask. Consult your physician if this occurs.
- If re-using the FlexiFit™ 432 Full Face Mask between different patients, refer to the Multi Patient Cleaning Instructions for the FlexiFit™ 432 Full Face Mask.
- Replace the Non-Rebreathing Valve every 12 months.
- Do not use the FlexiFit™ 432 Full Face Mask Elbow on other Fisher & Paykel masks.
- Use the FlexiFit™ 432 Full Face Mask only with the headgear provided.

CAUTION

- Use the mask for its intended use only as directed in this booklet.
- US Federal Law restricts this device to sale by, or on order of a physician.
- Operating temperature range is 40-120°F (5-50°C).

CONTRAINDICATIONS

- This mask should not be used by anyone incapable of removing the mask unaided, without continuous supervision.

US Patent Numbers 6615834 / 6789541 / 6701926 / 6951218
Other Patents Pending

BIAS FLOW INFORMATION


The FlexiFit™ 432 Full Face Mask has exhalation ports to expel the air that you exhale from the mask.

It is important that these exhalation ports are not blocked by any object. This controlled leak ensures all exhaled CO₂ is expelled from the mask.

Pressure (cm H ₂ O)	3	5	7	9	11	13	15	17	19	21	23	25
Bias Flow (l/min)	20	25	28	33	37	40	43	46	49	52	55	57

RESISTANCE TO FLOW

Pressure drop through mask @ 50 l/min 0.3 cm H₂O

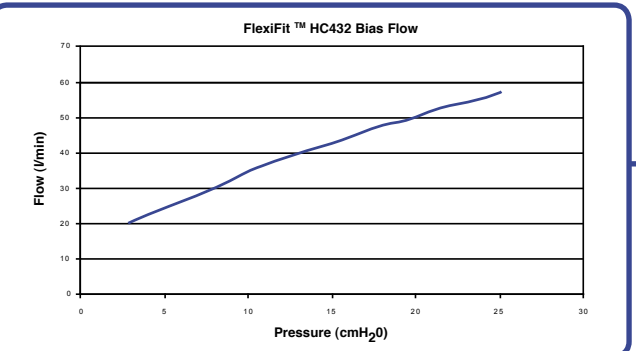FlexiFit™ 432 Full Face Mask Multi-Patient Cleaning Instructions

Your FlexiFit™ 432 Full Face Mask is provided patient ready. To reuse between patients please follow the guidelines below.

CLEANING / DESOILING

Before disinfecting, disassemble all components from the mask base. Clean all surfaces of mask and components thoroughly with a mild detergent solution and a soft bristled brush. Ensure all visible soil is removed. Rinse thoroughly under running potable water to remove detergent solution.

FAILURE TO PROPERLY CLEAN THIS DEVICE MAY RESULT IN INADEQUATE DISINFECTION

HIGH LEVEL DISINFECTION PARAMETERS

MASK FRAME	THERMAL DISINFECTION
	Soak parts in hot water bath for: 30 minutes at 70°C (158°F) or 30 minutes at 75°C (167°F)
	OR CIDEX® OPA Follow manufacturer instructions for high level disinfection. Immerse for 30 minutes at 20°C (68°F). Rinse thoroughly in potable water.
	OR AUTOCLAVING Gravity displacement type autoclave. Components should be enclosed in a sterilization peel pouch for 30 minutes exposure at 121°C (250°F).
	OR STERRAD® 100S Follow manufacturer instructions for Sterrad® 100S Sterilization system.
	WASH IN WARM SOAPY WATER AND LEAVE TO DRY
	REPLACE WITH NEW PART

The FlexiFit™ 432 Full Face Mask has been validated to withstand 20 cycles of the above parameters. Do not use if badly damaged.

Disclaimers: It is the user's responsibility to qualify any deviations from the recommended methods of reprocessing. Some discoloration may occur after re-processing.

Contact your local healthcare dealer for additional information.

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